trophon[®] Printer

User Manual

trophon[®] Printer

User Manual

Read this manual before operating the trophon® Printer to ensure correct procedures are followed.

For further information, contact your customer service representative or visit the Nanosonics website.

©2024 Nanosonics Limited. All rights reserved.

The content of this manual are correct at the time of product sale.

Nanosonics' technologies are protected by a combination of patents, trademarks, and exclusive proprietary rights.

trophon® is a registered trademark of Nanosonics Limited.

Conforming to Low Voltage Electrical Equipment Directive (EU) 2014/35/EU & Electromagnetic Compatibility Directive 2014/30/EU

Conforming to UK Electrical Equipment Safety Regulations 2016 (SI 2016/1101) and UK Electromagnetic Compatibility Regulations 2016 (SI 2016/1091)

www.nanosonics.com

Table of Contents

1.	INTRODUCTION AND SETUP	4
2.	INSTRUCTIONS FOR USE	6
3.	ROUTINE CARE AND MAINTENANCE	6
4.	TROUBLESHOOTING	7
5.	PRINTER DISPOSAL	9
6.	PRINTER TECHNICAL SPECIFICATIONS	9
APPENDIX A.....		100

1. INTRODUCTION AND SETUP

The trophon Printer provides a label record of the most recent high level disinfection cycle directly from the trophon device. The label is generated at the completion of each high level disinfection cycle and can be placed directly on to a patient medical record or other documentation. Printed label records remain legible for up to 7 years if stored as per the specifications outlined on the label roll packaging.

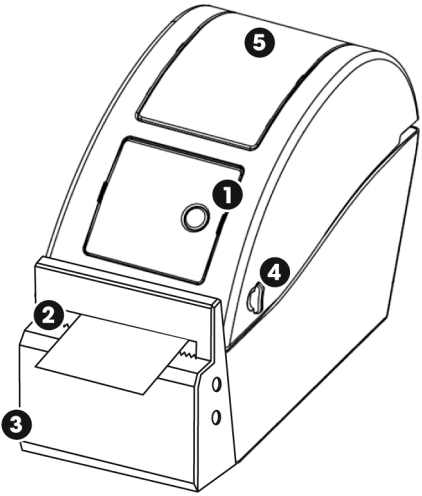
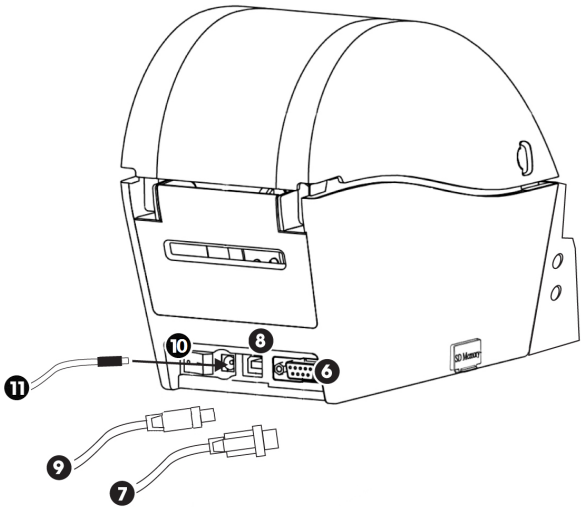
The printer is compatible with the trophon EPR and trophon2 devices. Accuracy and life time of records is not guaranteed if an unapproved accessory (printer or label roll) is used or if records are stored outside the storage specifications.

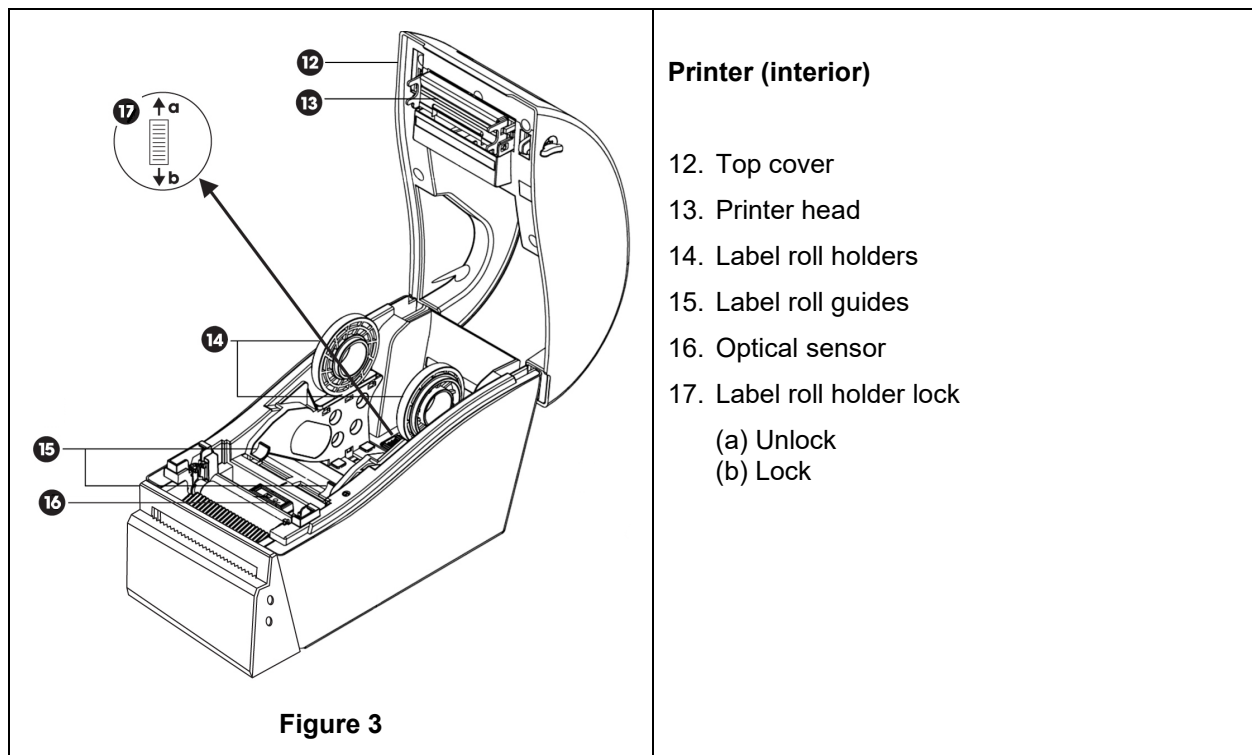
For any operation, fault or maintenance queries, please contact your customer service representative.

1.1 Important Labels, Symbols and Warning

	Caution		Consult instructions for use
	Fragile / Handle With Care		This Way Up
	Legal Manufacturer		Date of Manufacture
	Separate collection for electrical and electronic equipment		Keep Dry
	Warning: Hot Surface		Warning: Moving parts
	Conforming to Low Voltage Electrical Equipment Directive (EU) 2014/35/EU & Electromagnetic Compatibility Directive 2014/30/EU		Restriction of Hazardous Substances
	Product Number		Serial Number
	EU Importer		European Authorised Representative
	Conforming to UK Electrical Equipment Safety Regulations 2016 (SI 2016/1101) and UK Electromagnetic Compatibility Regulations 2016 (SI 2016/1091)		

1.2 Overview of Printer Features

 Figure 1	Printer (front) 1. Feed button 2. Paper exit channel 3. Auto cutter 4. Top cover open lever 5. Media view window
 Figure 2	Printer (back) 6. Serial cable interface 7. Serial cable 8. USB cable interface 9. USB cable 10. Power switch 11. Power cable

1.3 Printer – Initial Set-Up

- Unpack the printer from its packaging
- Position the printer securely on a bench beside the trophon device. The trophon Printer Wall Mount or trophon Printer Cart Mount accessories should be used to mount the printer to the wall or cart.
- Follow the setup instructions in Appendix A for powering on the printer and connecting it to your trophon device.

NOTE: Ensure the printer power cable is connected to the printer before connecting to an earthed mains outlet or switching power switch to on position.

NOTE: Keep the label roll holders in the “unlock” position. Please note the orientation of the label roll holders in Figure 3, as they can be rotated. Incorrect orientation will prevent correct loading of the label roll into the printer.

2. INSTRUCTIONS FOR USE

- Once the printer is connected to the trophon device via the appropriate cable provided (either USB or serial cable, depending on the port available on your trophon device) and powered on, the printer will be automatically detected.
- It is recommended to complete a test print to confirm successful connection to the device.
- Test prints and all other printer settings can be accessed on the trophon device screen by selecting either *Menu* → *Printer* or *System Setup* → *Printer Settings*. Please check your trophon device.

3. ROUTINE CARE AND MAINTENANCE

3.1 Replacing printer Label Roll

The trophon device will display an onscreen notification when the label roll requires replacement.

Follow the setup instructions included in Appendix A to load the label roll.

3.2 Cleaning the Printer

- Turn off the printer and unplug from the power socket before cleaning.
- Do NOT submerge the printer, or pour liquids over it.
- Keep the power socket completely dry (see Figure 2).
- Remove the label roll paper from the printer prior to cleaning.
- Clean the printer feeder area routinely to ensure there is no glue or paper residue. Use tweezers to remove any excess residue from the paper cutter area.
- For cleaning, wipe the outside covers of the printer with a cloth moistened with a mild, soapy solution, until all surfaces are visibly clean, taking care to ensure that liquid does not come into contact with the power socket, printer interior, or label roll.
- Clean the printer head when replacing the label roll using a soft cloth with ethanol. Refer to Figure 3 for the printer head location.
- For disinfecting, wipe all external surfaces of the printer with an Isopropanol or Quat Wipe.

4. TROUBLESHOOTING

4.1 Connecting and Disconnecting the Printer

The printer is automatically detected and enabled when connected to the trophon device.

The printer can be enabled or disabled from the printer settings accessed on the trophon device screen.

4.2 Printer – Error Messages

Printer error messages are displayed on the screen of the trophon device. If an error message is displayed, two options are available:

1. *Retry*

- After the error has been addressed, select *Retry*
- The trophon will attempt to re-print the last disinfection record.
- If the last disinfection record did not reprint, the two options *Retry* and *Ignore* will again display.

2. *Ignore*

- If the printed record is NOT required and you do NOT want to *Retry* the print, select *Ignore*.

NOTE: If *Ignore* is selected, the disinfection label record will not be printed.

4.3 Printer - Errors

Error	Check for the following
No Power to the printer	The printer is plugged in to power and switched on. The printer displays a small green light when powered on.
Printer LED Indicator – solid Red light	The paper cover is closed securely. Printing error - printer head open or cutter error.
Printer LED Indicator – flashing Red light	Printing error - head open, paper empty, or paper jam.
Printer LED Indicator – Amber light	The system is clearing data from printer.
Printer LED Indicator – flashing Green light	Printer is paused. Press the button on the Printer to re-start the Printer.
Printer Not Detected	The printer is plugged in correctly to the trophon device. The power is switched ON at the wall outlet and the printer. The printer is a Nanosonics approved printer - a trophon logo is printed on the front screen
Printer Not Compatible	The printer is a Nanosonics approved printer - a trophon logo is printed on the front screen
Printer Paused	Select <i>Retry</i> on the trophon device screen.
Printer Out of Paper	Follow the setup instructions in Appendix A for loading a new label roll.
Printer Paper Jam	Labels may be stuck inside the printer mechanism. Check and remove any jammed labels. Once cleared, select <i>Retry</i> on the trophon device screen.
Printer Cover Open	Close the printer cover. Select <i>Retry</i> on the trophon device screen.
Printer Cutter Error	Labels may be stuck inside the printer. Check and remove any jammed labels. Once cleared, select <i>Retry</i> on the screen.
Printer Head Overheat	Select <i>Retry</i> on the screen. If the error recurs, switch the printer off and allow cooling, then switch the printer on again and select <i>Retry</i> on the screen. Ensure the printer is positioned in a ventilated area and not confined space.
Printer Error	Switch the printer off and on and select <i>Retry</i> on the screen.
No Print on Label	Check the label roll is loaded correctly and the printer is switched on. Hold the feed button down to ensure the label roll is fed through properly. A test print is recommended when loading a new label roll into the printer to ensure correct loading.

5. PRINTER DISPOSAL

The following information is for EU-member states only:

This symbol indicates that the product must not be treated as general waste and must be disposed of correctly. For more detailed information about the collection and recycling of this equipment, please contact your supplier.

6. PRINTER TECHNICAL SPECIFICATIONS

Electrical Specification – printer	Input: DC 24V, 2A
Electrical Specification – External universal switching power supply	Input Voltage: 100 – 240Vac Input Current: 1.8A, 50/60Hz Mains Inlet: IEC type C13 Equipment must be connected to an earthed outlet using power cable supplied. Output: DC 24V, 2.08A
Environmental Specification	Operation: 41-104°F (5-40°C), 25 – 85% non-condensing Storage: -40 – +140°F (-40 – +60°C), 10 – 90% non-condensing
Physical Dimension	9.84" (L) x 4.33" (W) x 6.89" (H) 250mm (L) x 110mm (W) x 175mm (H)
Data Interface	RS-232C serial cable (max 115,200 bps)
Electromagnetic Compliance	FCC Class B, CE Class B, C-Tick Class B

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Nanosonics Europe GmbH
Nanosonics Limited UK