

trophon²

Protection by design



High-level disinfection for ultrasound probes

Consumables and accessories



With AcuTrace[®] (RFID) technology to enhance digital record-keeping and deliver easy, reliable workflows for ultrasound probe high-level disinfection (HLD).

nanosonics
Infection Prevention. For Life.



High-level disinfection device

trophon®2

The latest innovation in ultrasound probe HLD featuring an enhanced design, streamlined workflows, plus AcuTrace RFID technology for supporting digital record traceability.



trophon NanoNebulant®

trophon NanoNebulant is a proprietary disinfectant liquid with 35% hydrogen peroxide chemistry to achieve effective high-level disinfection of the entire ultrasound probe including the shaft and handle. Incorporates AcuTrace technology to store batch and expiry details directly onto your trophon2 device.



trophon Chemical Indicator


The trophon Chemical Indicator responds to a combination of specified time, disinfectant dosage and disinfectant concentration. These critical parameters need to be met for the chemical indicator to change to a pass colour.

This qualitative colour change provides independent confirmation of success of each disinfection cycle. Indicator packaging now incorporates AcuTrace technology to record and store batch and expiry details directly into trophon2.



trophon AcuTrace Operator Card

Incorporates AcuTrace technology to digitally capture operator interaction details, stored directly on your trophon2.

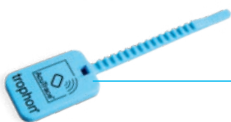

trophon AcuTrace Medical Instrument Tag

Incorporates AcuTrace technology to digitally capture ultrasound probe details, stored directly on your trophon2.



trophon Printer

The trophon printer is a convenient traceability solution, printing high quality trophon disinfection labels each and every cycle.



trophon Printer Label Roll

High quality printer labels with a minimum 7-year life span, appropriate for medical records.



trophon Logbook

Simplifies and systematises record-keeping while providing evidence of quality HLD reprocessing at health care facilities. Compatible with trophon Printer Labels.

Point-of-care set-up accessories



trophon Wall Mount 2

A convenient option for installing trophon2 where there are space constraints.


trophon Printer Wall Mount[^]

Custom designed for secure, horizontal mounting of the trophon Printer to a wall. Includes a cable tidy area for convenience.



trophon Cart

Makes trophon fully mobile for convenient point of care use and adaptability where there are space constraints.

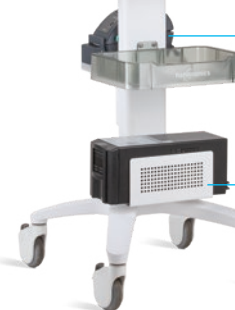

trophon Printer Cart Mount 2[^]

Designed for securely attaching the trophon Printer to the trophon Cart, in an easily accessible, convenient location.



trophon Wireless Ultrasound Probe Holder^{**}

Designed as a reusable accessory to facilitate high-level disinfection of wireless ultrasound probes. Integrates into your existing trophon device workflow and does not require any assembly, installation or hardware modification to your device.

Drying wipes



trophon Companion Drying Wipes

Used to dry probes before and after high-level disinfection. trophon Companion Drying Wipes are highly absorbent to leave surfaces completely dry.

trophon Companion Wipes Wall Mount

For wall-mounting trophon Companion Drying Wipes for easier access and dispensing.

Ultrasound probe storage cover



Clean Ultrasound Probe Cover

Custom designed probe covers to prevent endocavitary and surface ultrasound probes from recontamination after high level disinfection and before they are used on the next patient.

[^]Printer not included.

^{*}UPS not included. Contact your local sales representative for compatible UPS specifications. Disinfection cycles should not be run while trophon device is unplugged from the mains power.

^{**}Currently for use with GE Healthcare Vscan Air TM only.



trophon2 System Product Descriptions

Product Name	Product Code	Quantity
trophon2	N05000-ANZ	1 Device

Accessories

trophon NanoNebulant	N05001	1 Box of 6 Bottles*
trophon Chemical Indicator	N05003	1 Box of 300 Chemical Indicators
trophon AcuTrace Operator Card	N05006	1 Box of 10 Cards
trophon AcuTrace Medical Instrument Tag	N05007	1 Box of 10 Tags
trophon Printer	N00048-ANZ	1 Printer
trophon Printer Wall Mount	N00105	1 Mount
trophon Printer Label Roll	N00049-ROW	1 Roll of 600 Labels**
trophon Logbook	N00098	1 Box of 5 Logbooks
trophon Wall Mount 2	N05009	1 Mount
trophon Cart	N05010	1 Box of 1 Cart
trophon Printer Cart Mount 2	N05013	1 Mount
trophon Cart UPS Bracket	N00110	1 Bracket
trophon Wireless Ultrasound Probe Holder (WUPH-A)	N05030	1 Unit
trophon Companion Drying Wipes	N00113	1 Box of 6 Cannisters x 100 Wipes
trophon Companion Wipes Wall Mount	N00114	1 Mount
Clean Ultrasound Probe Covers	N00102	1 Box of 100 Covers

*Minimum order quantity: 2 boxes.

**Parametric Release Label: 1 Roll of 350 labels.

For medical and health professionals only

Nanosonics Limited
(Manufacturer)

7-11 Talavera Road
Macquarie Park, NSW 2113, Australia
T: +61 2 8063 1600
E: info@nanosonics.com.au
W: www.nanosonics.com.au



Protection by design



Infection Prevention. For Life.